


Lecture 25 – Other Imaging modalities


This lecture will cover:

- Optical Imaging
- Near-infrared Spectroscopy (fNIRS)

How important is optical imaging?

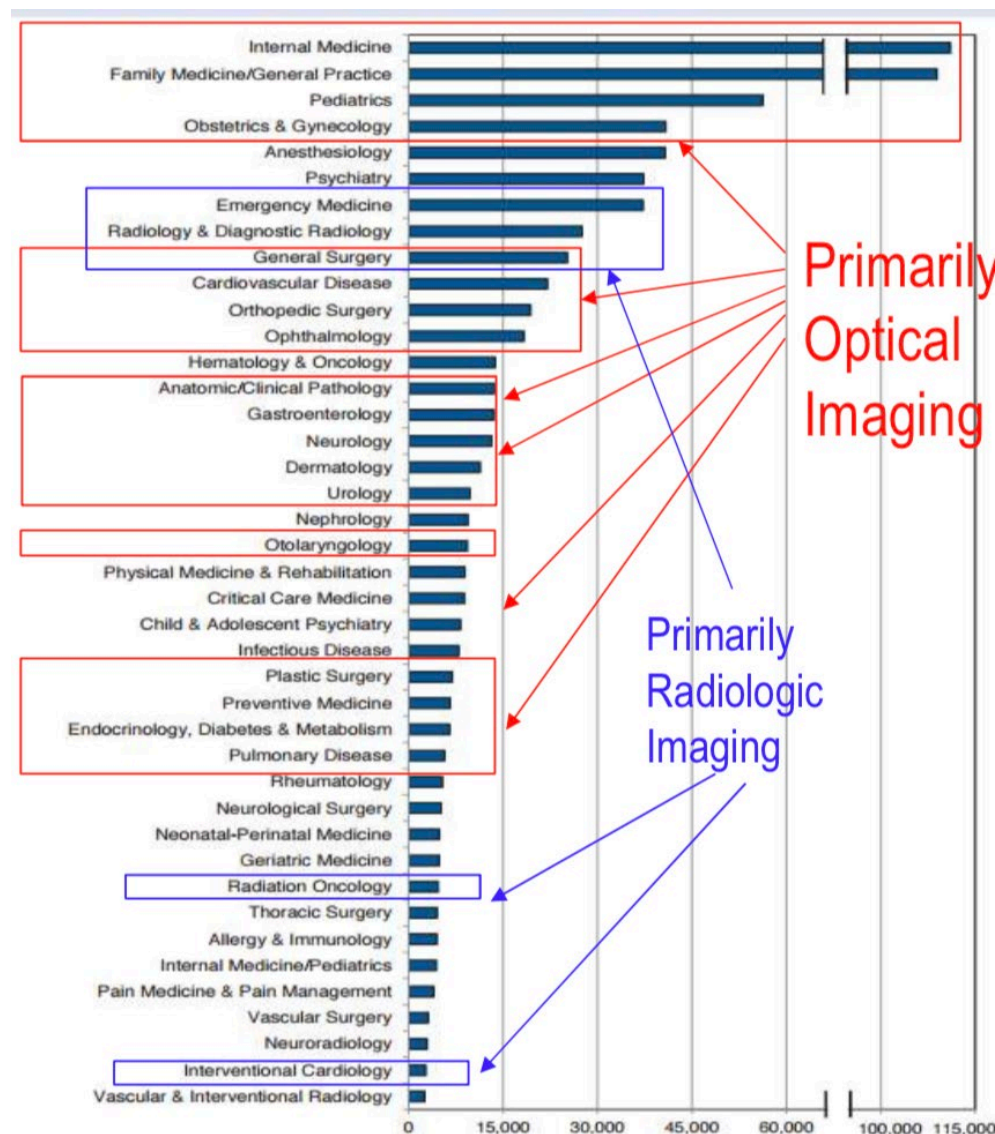


Number of Active doctors

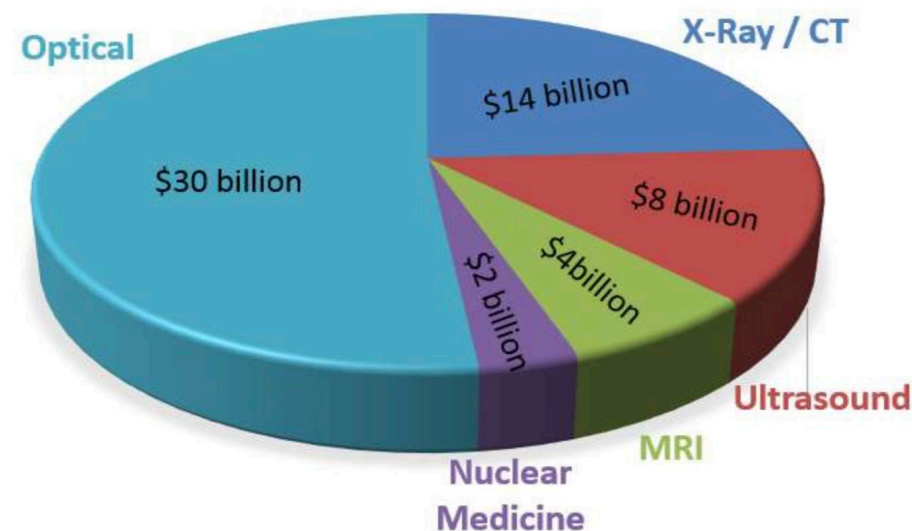

Market analysis

Refs: BCC Research, IBIS World, AMA



Optical imaging



上海科技大学
ShanghaiTech University


Advantages:

High sensitivity

Rich contrast

Non-ionizing radiation

Low cost

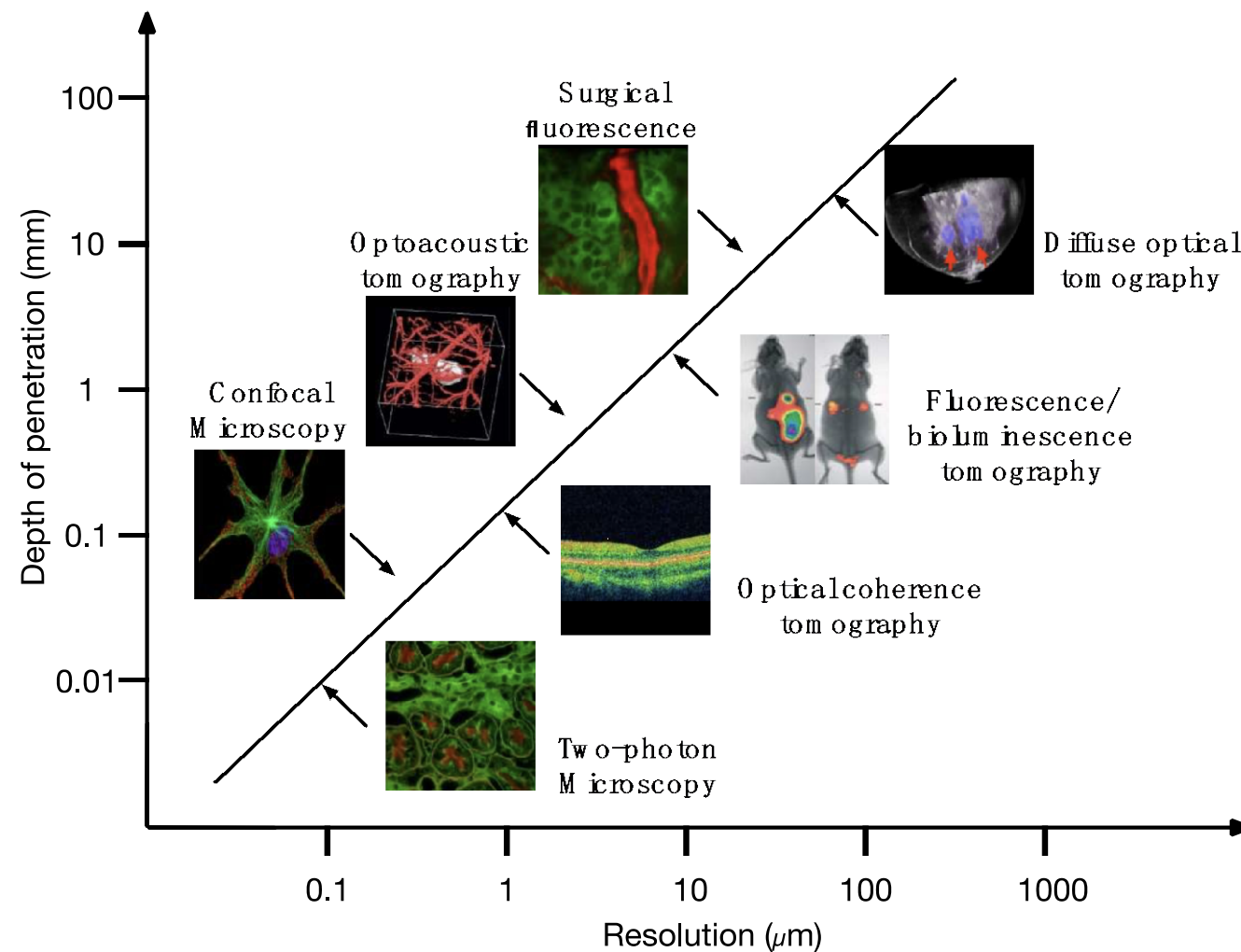

立志成才 报 国 裕 民




Confocal Microscopy (共聚焦显微镜)



上海科技大学
ShanghaiTech University

Confocal microscopy

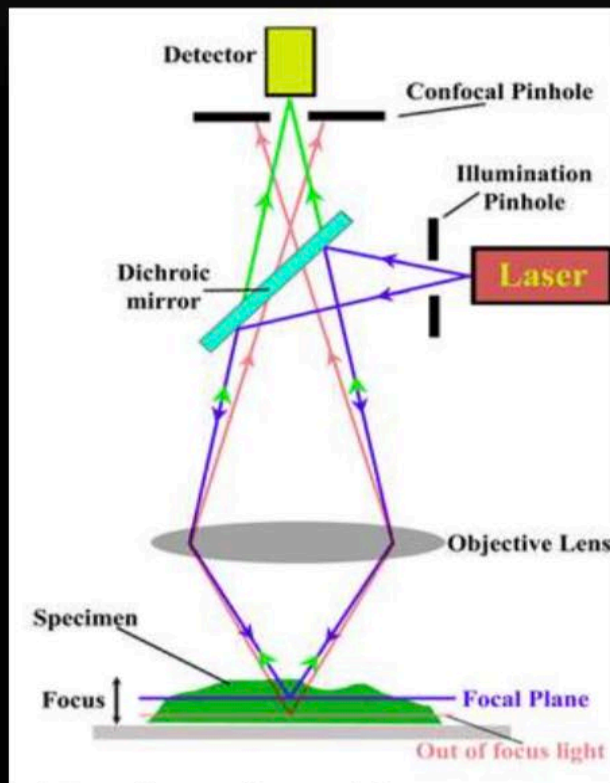

<https://www.jic.ac.uk/>

In vitro




microscopy.nri.ucsb.edu

Resolution $\cong 0.20 \mu\text{m}$
Penetration $\cong 40 \mu\text{m}$

In vivo



<http://www.confocaltraining.com>

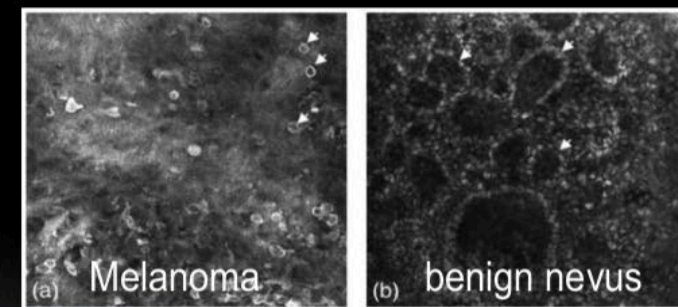



Resolution $\cong 2.0 \mu\text{m}$
Penetration $\cong 400 \mu\text{m}$

■ Super resolution Microscopy (超分辨显微镜)



上海科技大学
ShanghaiTech University



- Stimulated Emission Depletion Microscopy (STED)
- Photoactivated localization microscopy (PALM), Stochastic optical reconstruction microscopy (STORM)

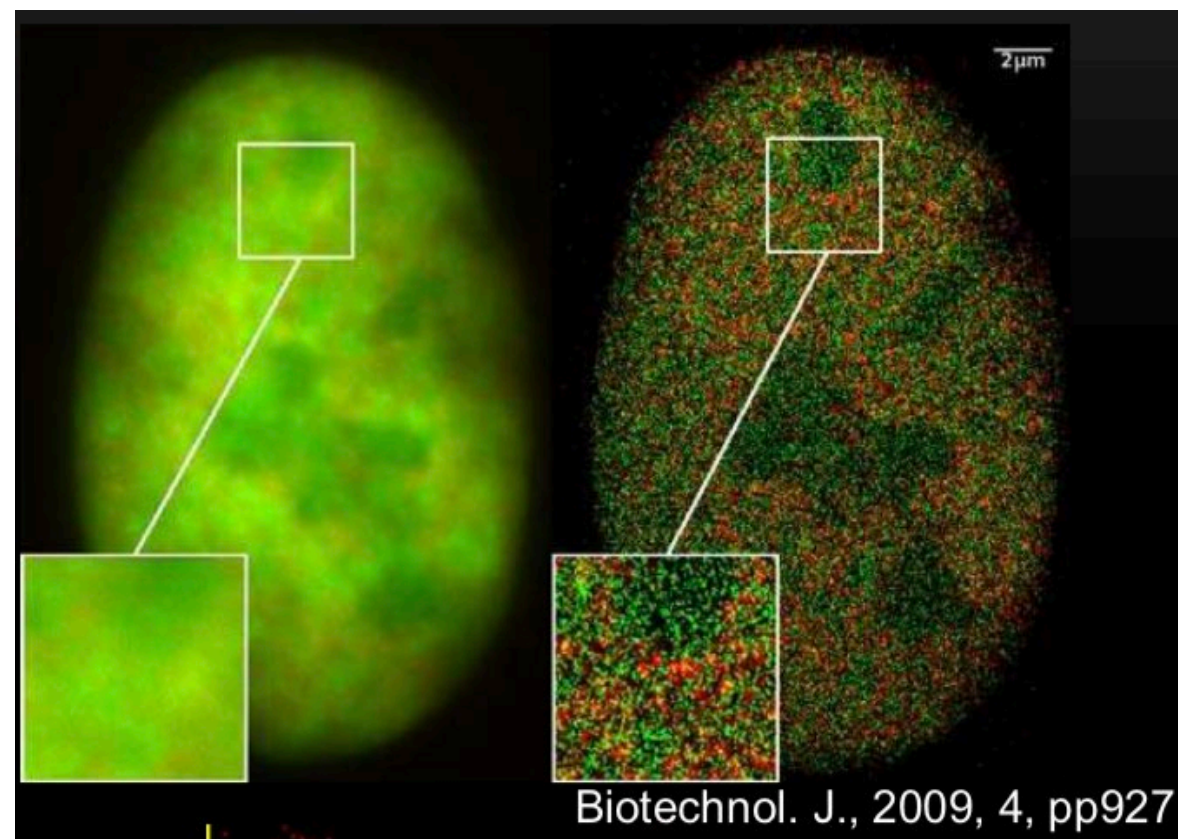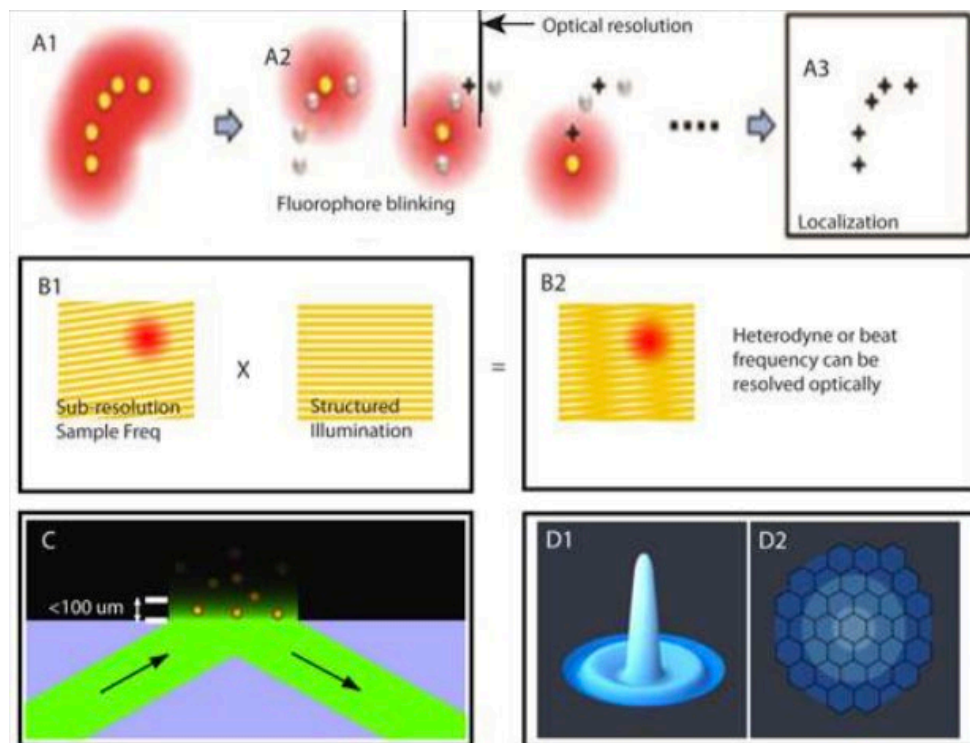


Optical Coherence Tomography (光学相干断层扫描技术)



上海科技大学
ShanghaiTech University

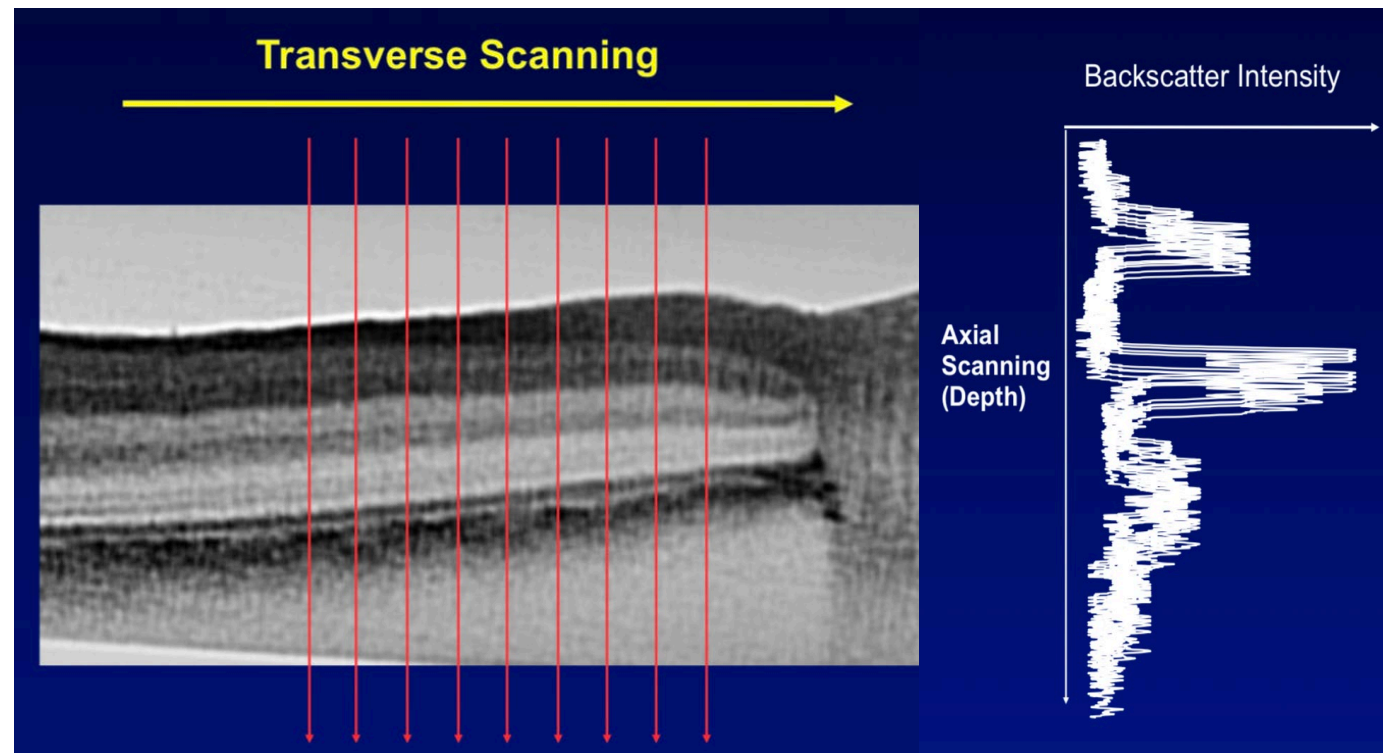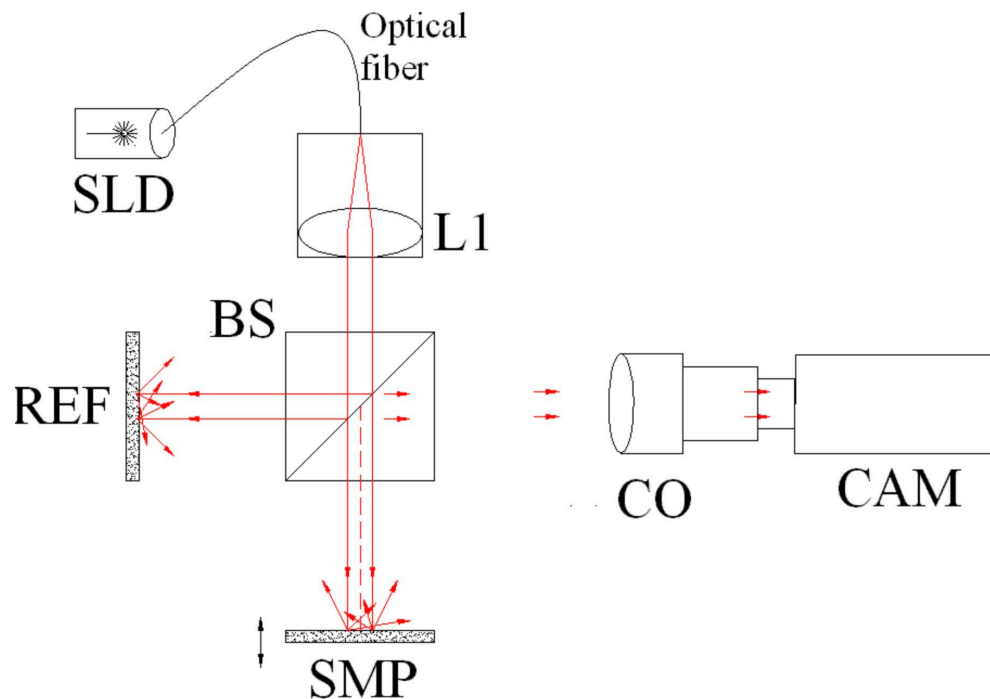

~1 um resolution (micron)


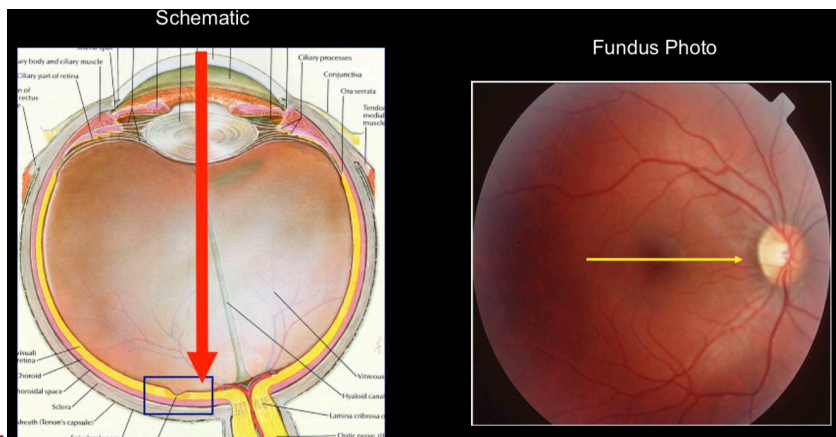


立志成才 报效国家

■ Endoscopy (内窥镜)



Resolution = ~ 0.10 mm
Penetration = ~ 2 mm

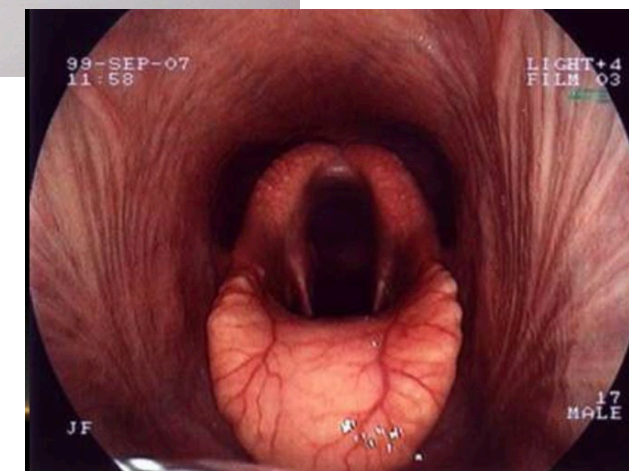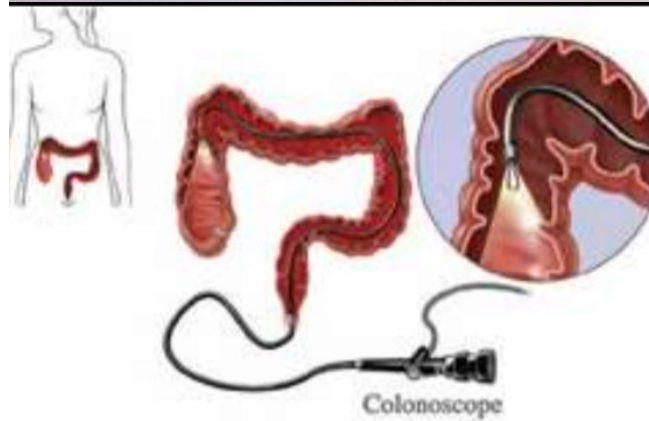

Introduction to fNIRS

- Near-infrared spectroscopy: irradiating the subject with light in the near-infrared region and examining the change in absorbance of light;
- Measuring the changes of cortical deoxygenated and oxygenated hemoglobin concentrations

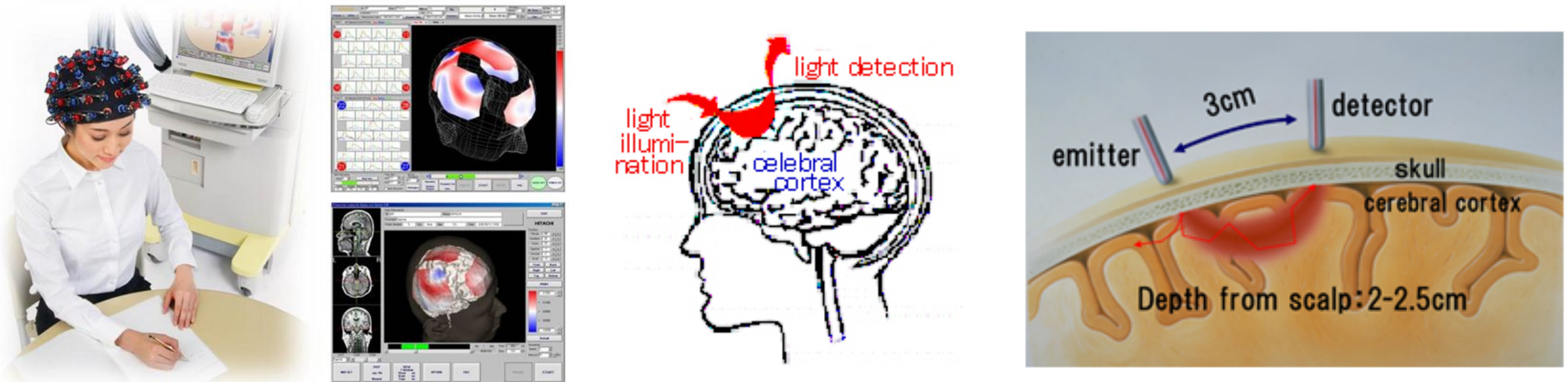


Fig. Demonstration of percutaneous measurement of cerebral cortex beneath the skull using fNIRS.

Optical Window

- 650–900 nm spectral interval;
- skin, tissue, and bone are mostly transparent;
- Hemoglobin (Hb) and deoxygenated-hemoglobin (deoxy-Hb) are strong absorbers

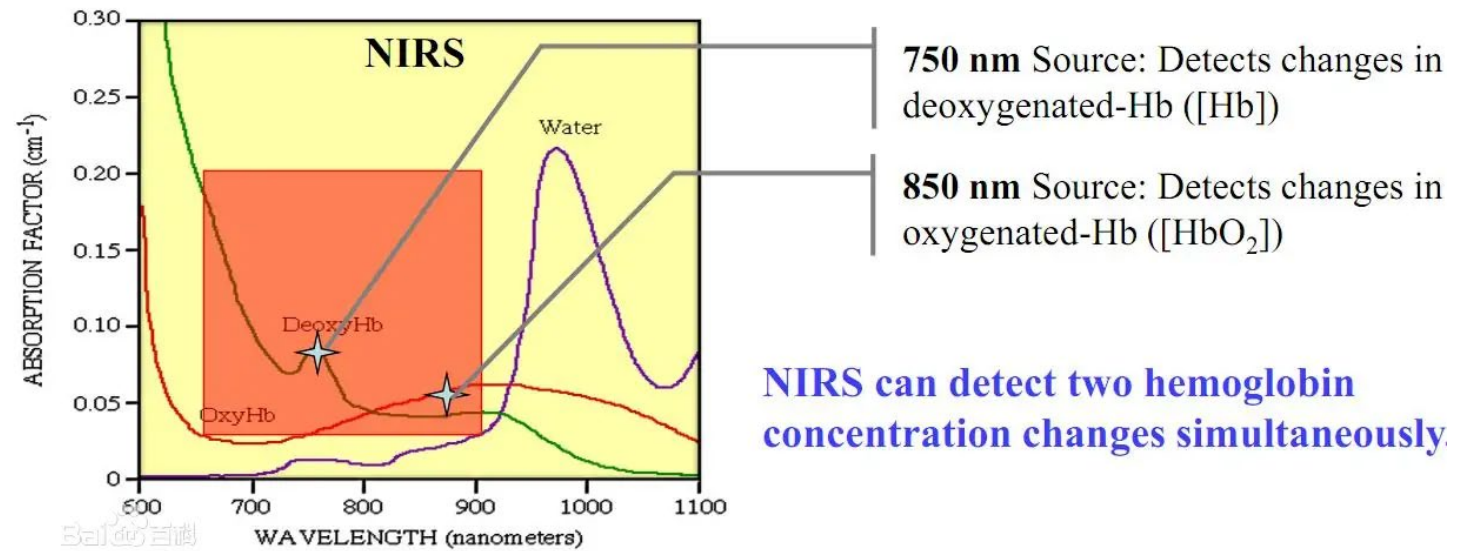
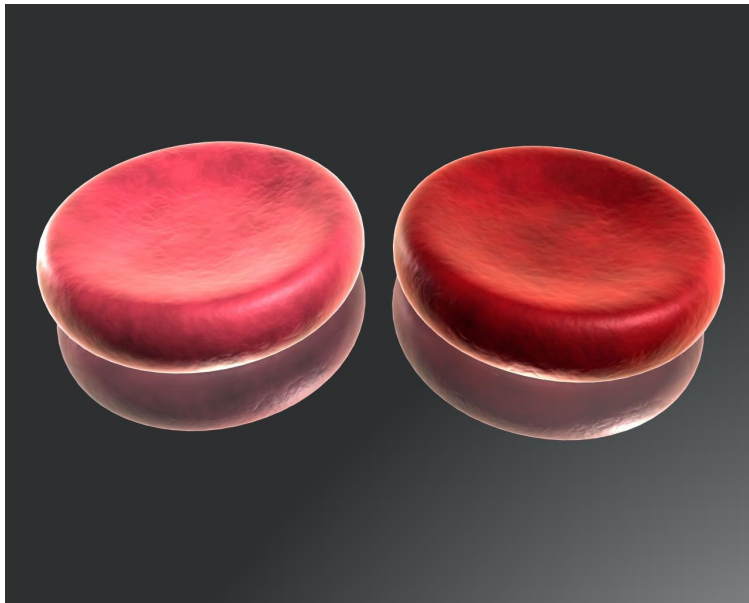


Fig. (Left) Oxygenated and Deoxygenated Hemoglobin; (Right) Absorption spectra for oxy-Hb and deoxy-Hb for near-infrared wavelengths.

Features

- Techniques: continuous Wave, frequency Domain, time-Domain
- Advantages: wearable and portable, noninvasiveness, low-cost modalities, perfect safety, high temporal resolution, full compatibility with other imaging modalities, and multiple hemodynamic biomarkers.
- Disadvantage: low brain sensitivity, low spatial resolution, and shallow penetration depth

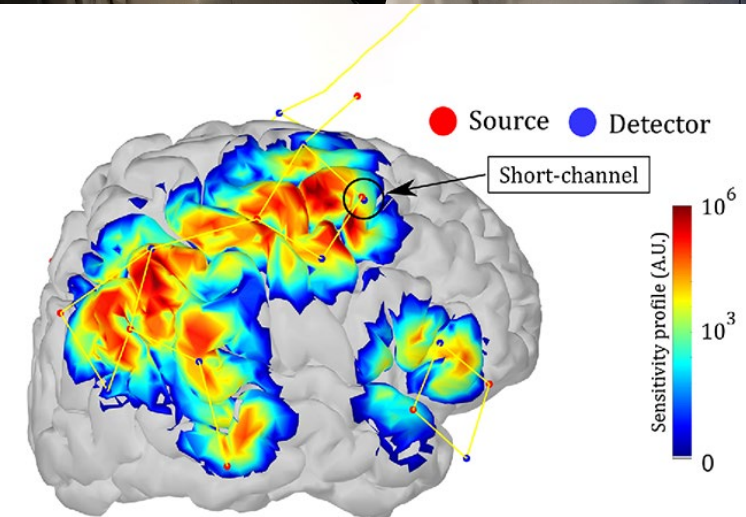
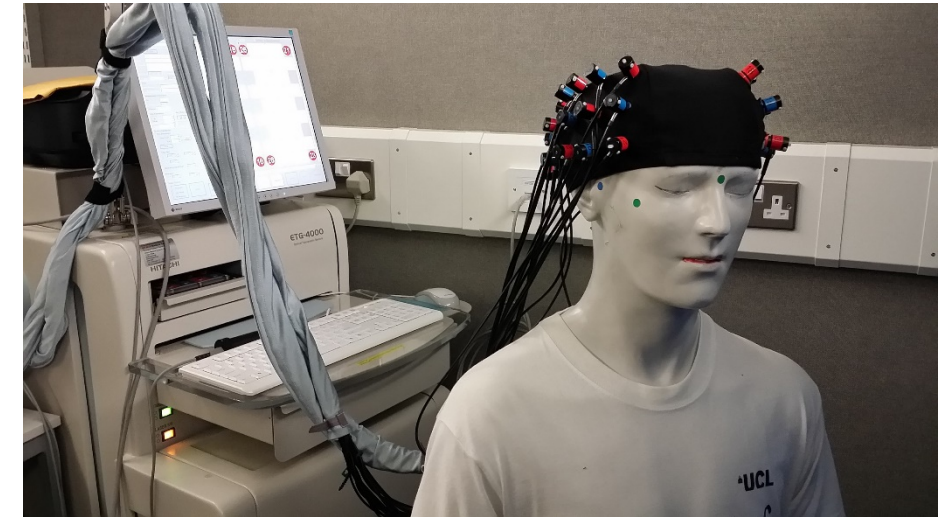


Fig. (Top) a fNIRS system. (bottom) fNIRS brain imaging